



ALLIANCE
BIO EXPERTISE



RIEGER

Industrievertretungen GmbH

www.rieger-iv.at

Wir sind für Sie da:

E-Mail: office@rieger-iv.at

Telefon: +43(1) 728 00 52


Fax: +43(1) 728 69 16

Rustenschacher Allee 10

1020 Wien, Austria

DOSYWEL[®]

Peristaltic pump

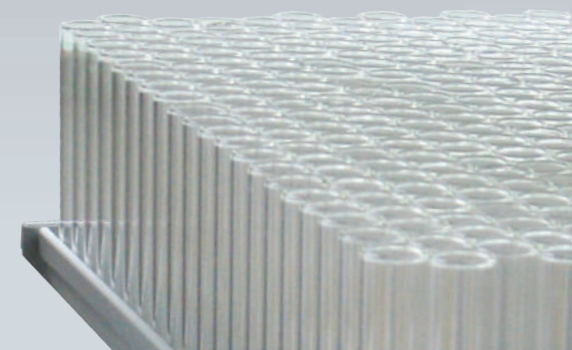



DOSYWEL®

Increase your productivity
quality
traceability

DOSYWEL®



- Up to 20 programs
- Editable:
 - ◆ Volume (1 to 9999ml)
 - ◆ Number of doses (1 to 999)
 - ◆ Time for interval (0.1 to 60.0 seconds)
 - ◆ Pump' speed (1 to 130% = 10 to 550 rpm)
 - ◆ Diluent (batch, name, volume)
 - ◆ Calibration mode (calibration starting the program, reminder each x sample).



Compact & light: 3.8kg
Easy transportable

QUICK & EASY CALIBRATION

Few steps and the calibration is done.
Calibration reminder: The DOSYWEL can also remind you every x sample, to perform a calibration

AVAILABLE OPTIONS :

- Footswitch



SPECIFICATIONS

DOSYWEL®



PERFORMANCES

Distribution volume	0.1 to 9999.9 ml
Pause time	1 – 99 seconds
Adjustable speed	1 – 130% 10 – 550 rpm
Connectivity	USB ports (Ticket printer, Footswitch, Keyboard, USB data export)
Communication	USB port. Ethernet
Internal tube dimensions	0.8 / 1.6 / 3.2 (standard) / 4.8 / 6.4 / 8.0 mm
Country of origin	100% France
Programs	20
Rotation way	Left or right
Foot switch	✓ (optional)

DIMENSIONS & ELECTRICAL SUPPLY

W x H x D	190 x 210 x 250 mm
Electrical supply	100-250v 50-60 Hz 10 A
Power supply	150 W

TRACEABILITY

Ticket printer	✓ (optional)
USB export	✓ (optional)

CONNECTED
MICROBIOLOGY



ALLIANCE BIO EXPERTISE



Empowered from our people's skills and experiences, we, ALLIANCE BIO EXPERTISE, together with our customers and partners decided to manufacture Microbiology automated equipment around essential values:
Customer Focus, Service, Loyalty, Innovations

**Ask your local partner
for a personal demonstration!**

ZI de courbouton , Le Tremplin – 35480 GUIPRY - FRANCE
E-mail : contact@abioexpertise.com
Phone : +33(0)2.40.51.79.53
alliance-bio-expertise.com